

# WHAT'S NEW?

## Tri-Ad™ 2.0

Adams Tricuspid Annuloplasty Band



### HOLDER



- End-to-end support with full edge accessibility for suturing
- High-visibility white with clearly marked implant size
- Visible semi-saddle aspect to accommodate the tricuspid annulus saddle
- Two prominent cut points that are easily visible and accessible

### SIZER



- Transparent with product pictogram for visualization of coverage and orientation
- Available in six sizes 26, 28, 30, 32, 34, 36 mm
- Added commissural recesses that are designed to avoid snagging of sutures
- Hole to facilitate use of a suture tag

### HANDLE



- Available in three lengths
- Securely attaches to sizer and holder
- Semi-malleable stem to accommodate varying anatomies



Pre-repair




Post-repair

Photos courtesy of Dr. David Adams, Mount Sinai Hospital, New York

**Medtronic**  
Further, Together

## Tri-Ad™ 2.0 Adams Tricuspid Annuloplasty Band, Model 900SFC

Order Number	Band Size (mm)	A (mm)	B (mm)	C (mm)	D (mm)	E (mm)	Orifice Area (mm <sup>2</sup> )
900SFC226	26	17.4	25.8	25.6	33.2	15.5	365
900SFC228	28	19.4	27.6	28.1	34.9	17.1	435
900SFC230	30	20.9	29.1	30.6	37.2	18.5	502
900SFC232	32	22.3	30.7	31.5	38.2	20.0	557
900SFC234	34	24.0	32.3	34.7	41.1	20.9	654
900SFC236	36	25.5	33.8	36.3	42.9	21.9	730

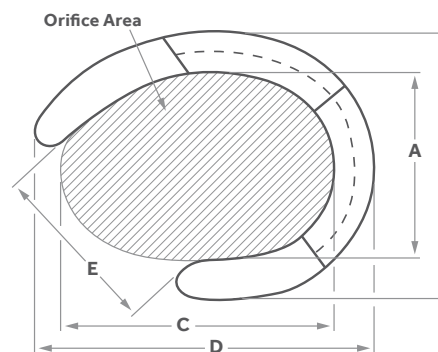

## Tri-Ad 2.0 Accessory and Ordering Information

Order Number	Components and Accessories
T7900P	Tri-Ad 2.0 sterilization tray
7900PS	Set of six Tri-Ad 2.0 sizers (26-36 mm)
7686 7686L 7686XL	Annuloplasty band and sizer handles

To order a full kit requires a set of sizers, a tray, and typically two of the shorter-length handles (the regular and long version will fit in the tray). The optional extra-long thoracoscopic handle will not fit into the tray and is sterilized separately.



For a listing of indications, contraindications, precautions, warnings, and potential adverse events, please refer to the Instructions for Use.

# Medtronic

Europe  
Medtronic International Trading Sàrl.  
Route du Molliau 31  
Case postale  
CH-1131 Tolochenaz  
Tel: +41 (0)21 802 70 00  
Fax: +41 (0)21 802 79 00

UC201901347EE © Medtronic 2018.  
All rights reserved. Printed in Europe.